

Reverse Preparation Technique


- Predictable reduction through the use of depth cuts
- Nearly perfect margin formation that is incredibly simple










Reverse Preparation Technique

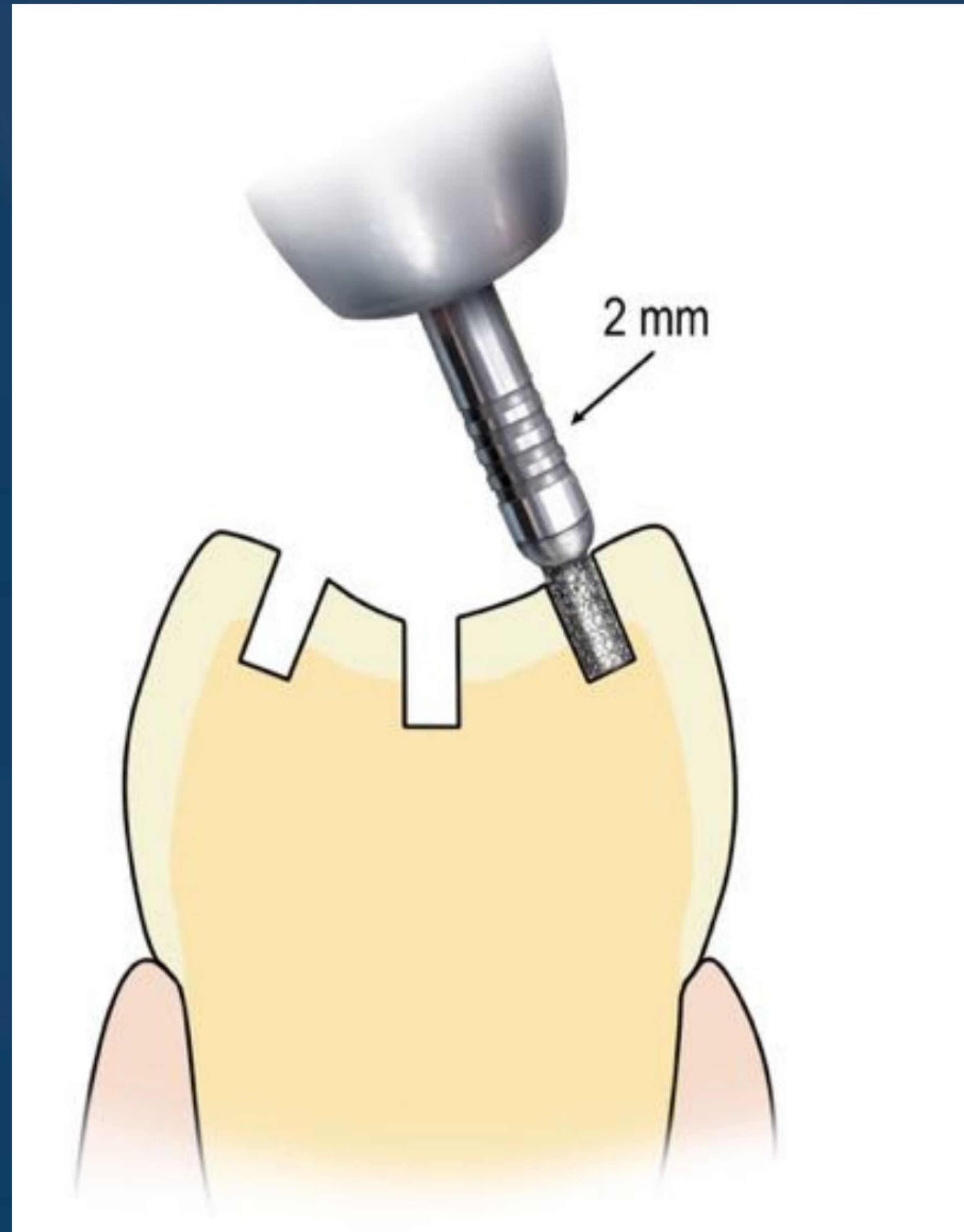


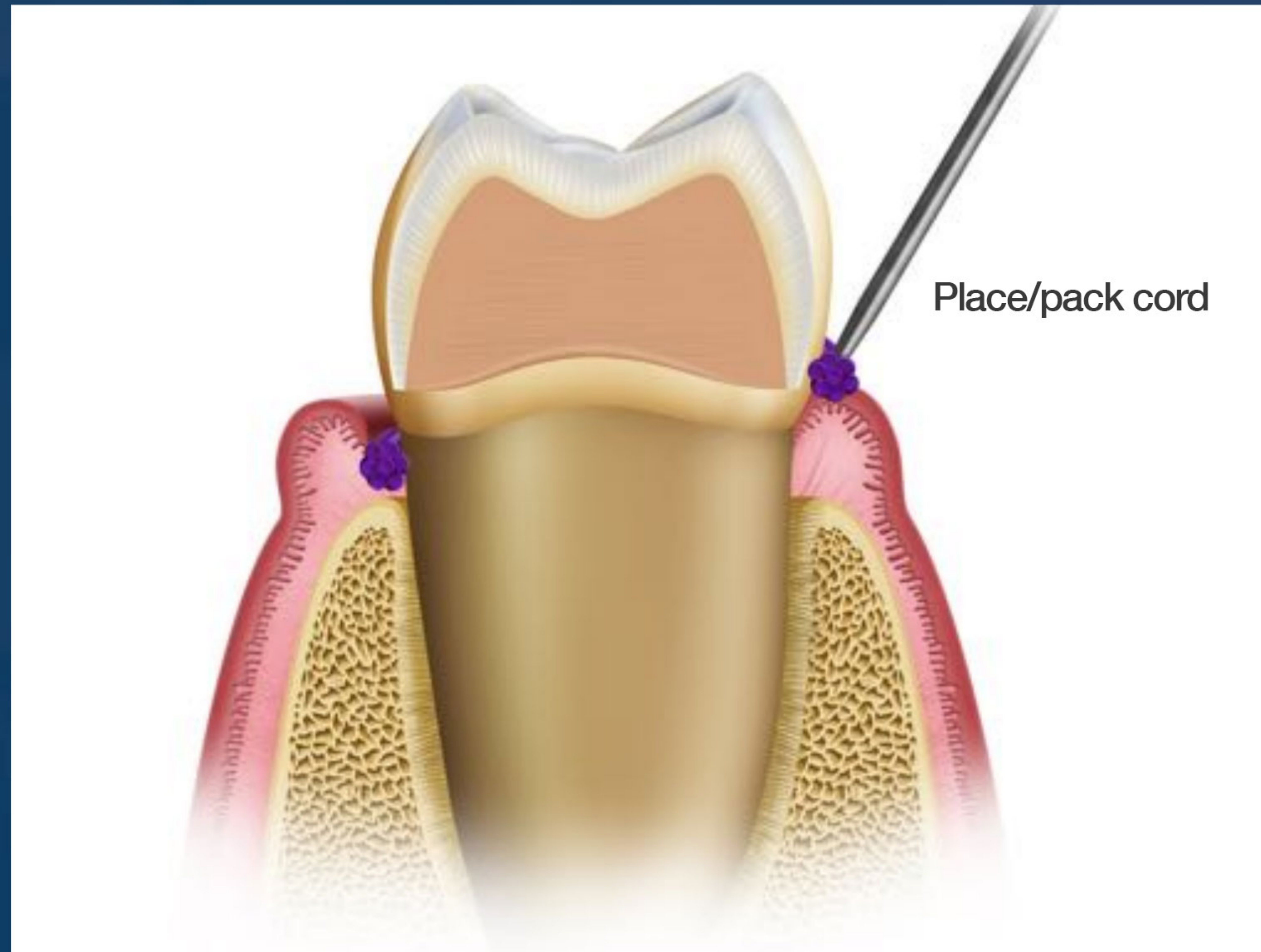





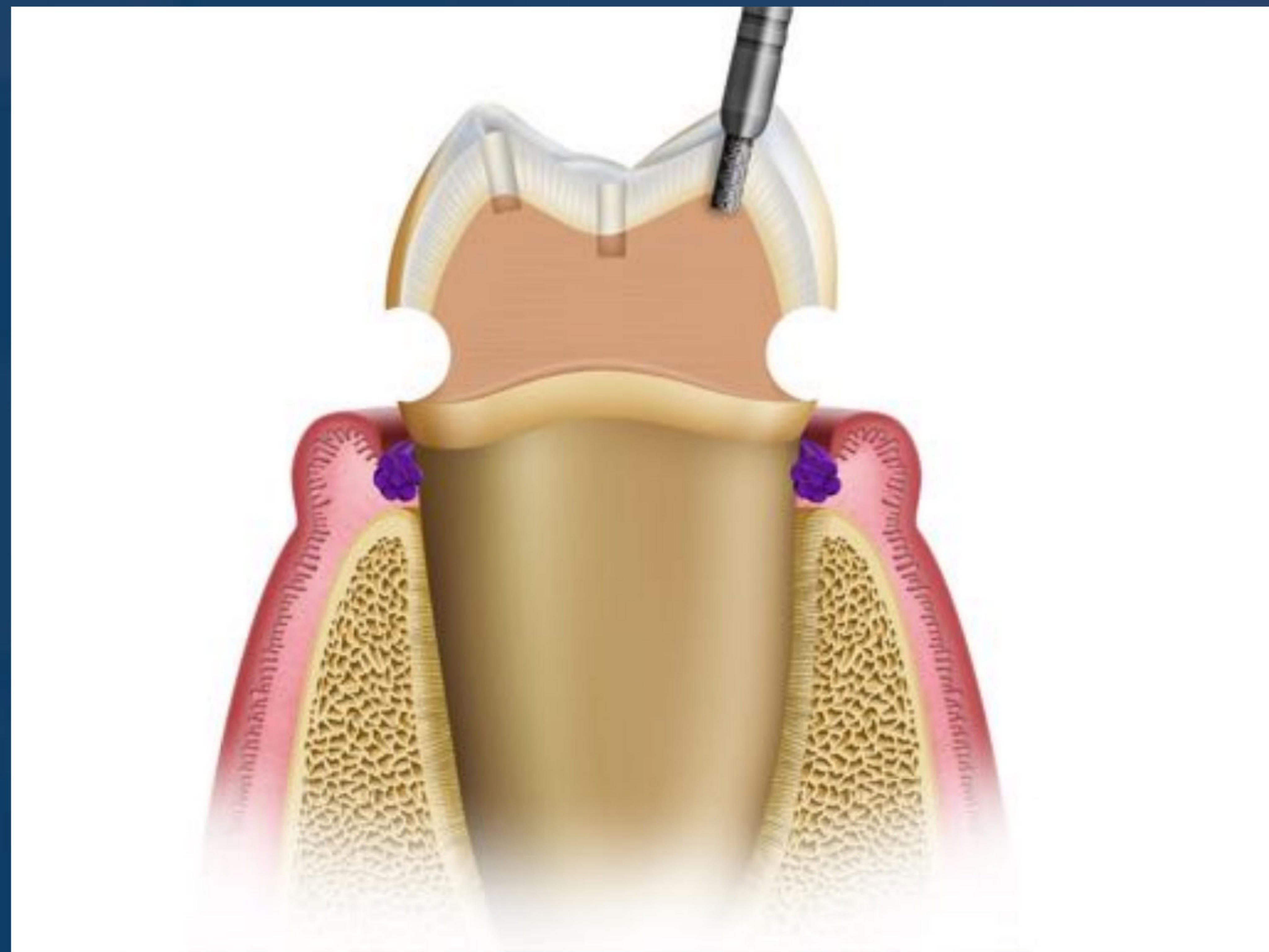


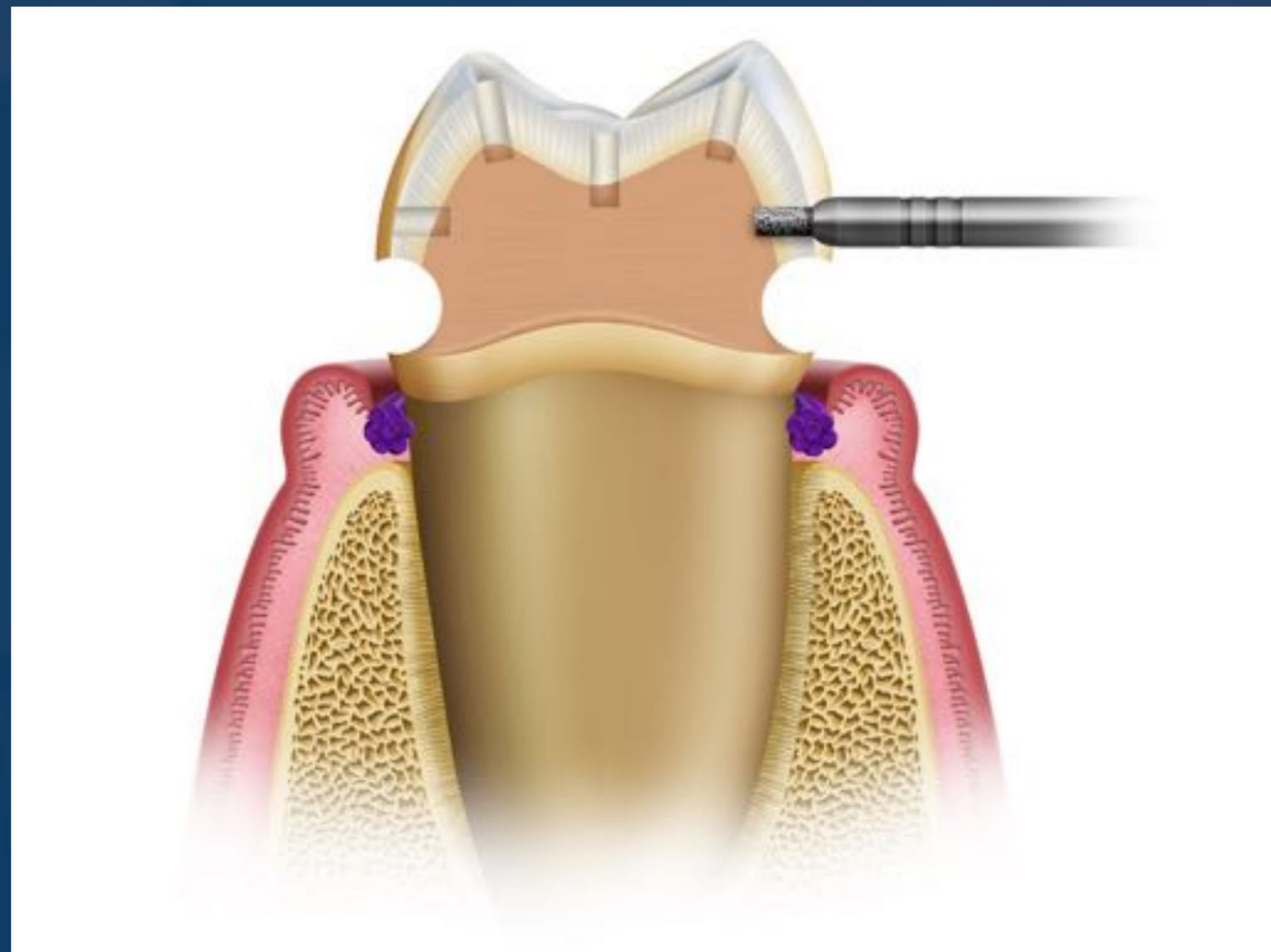


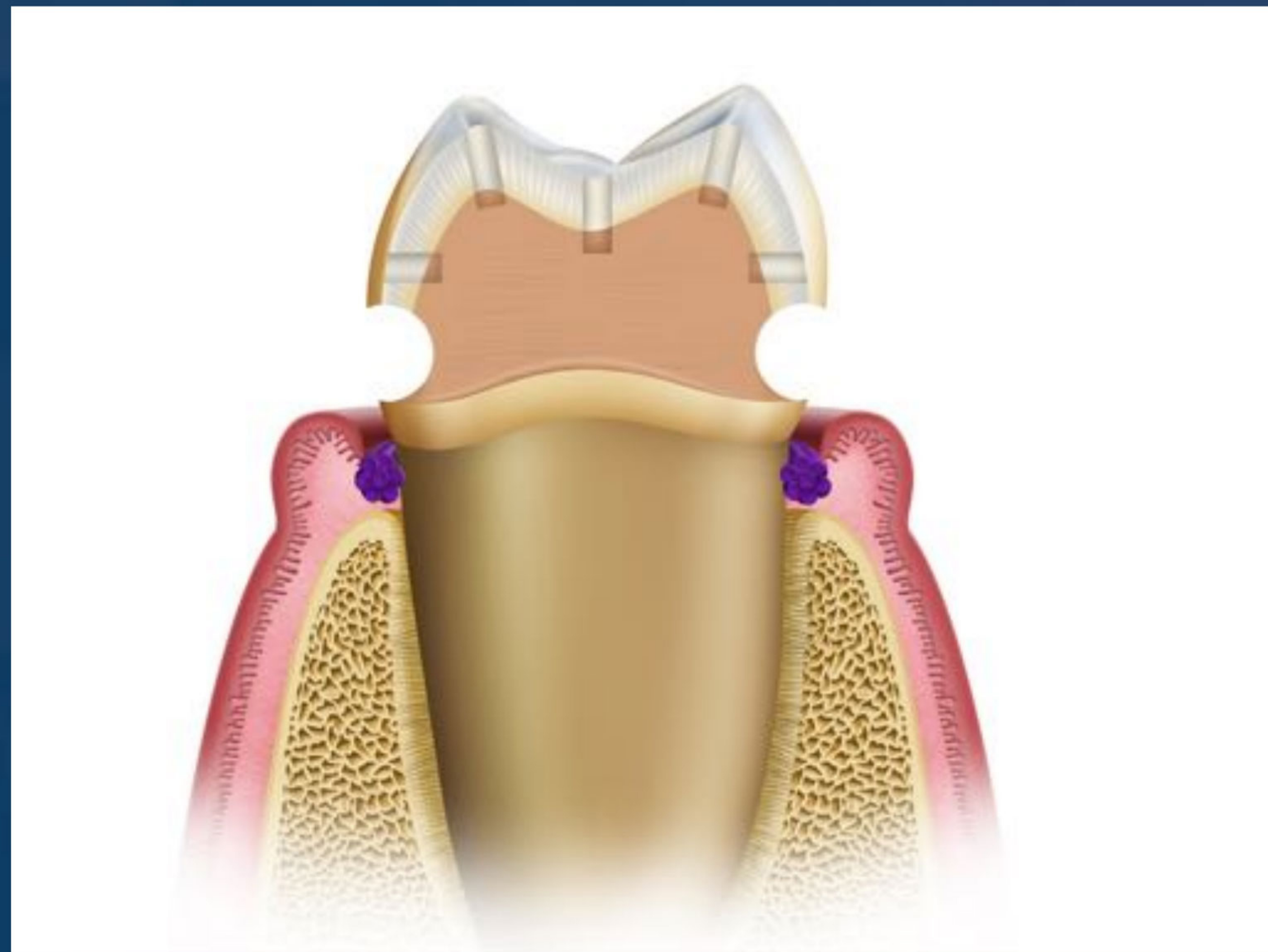


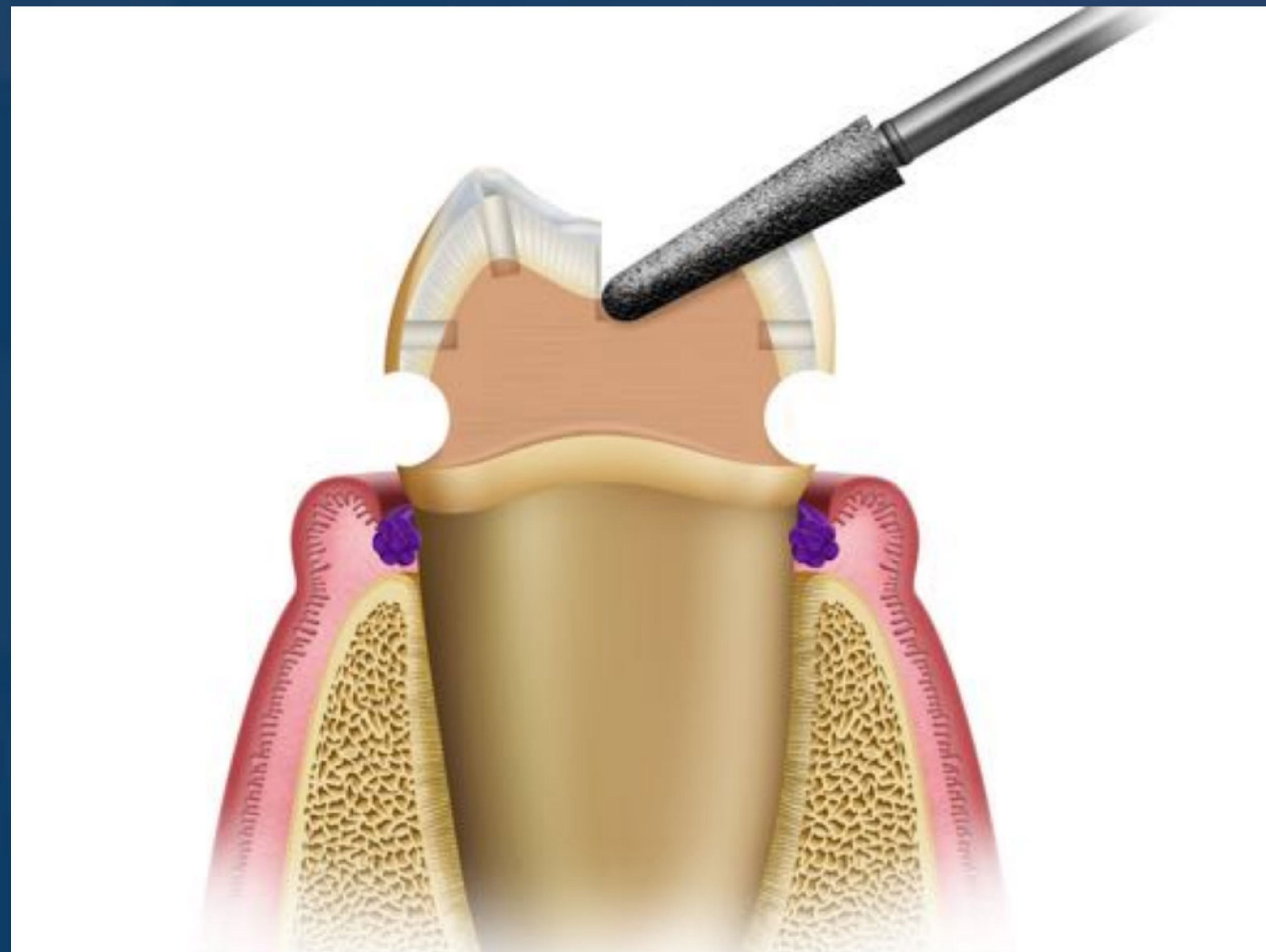


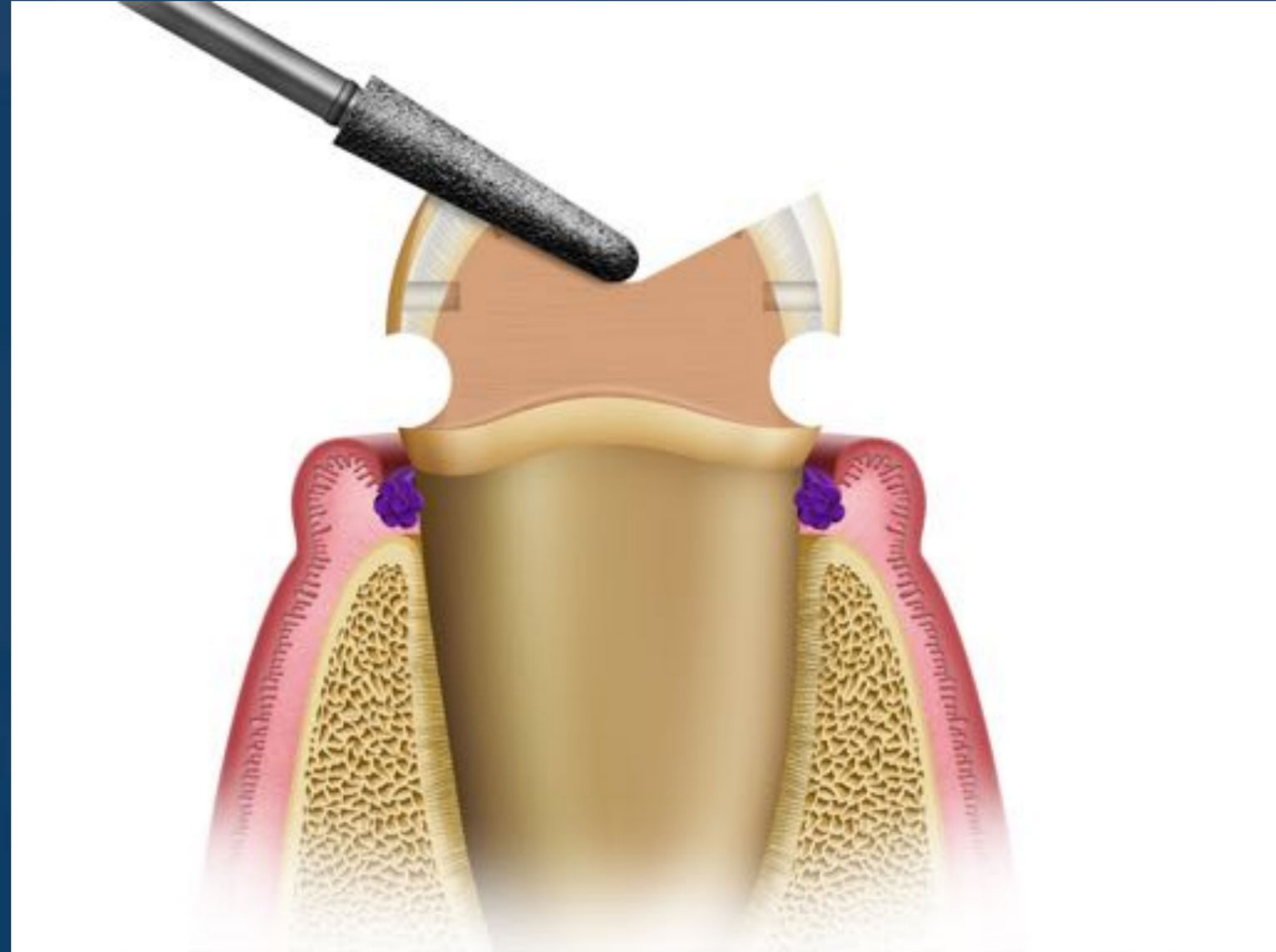

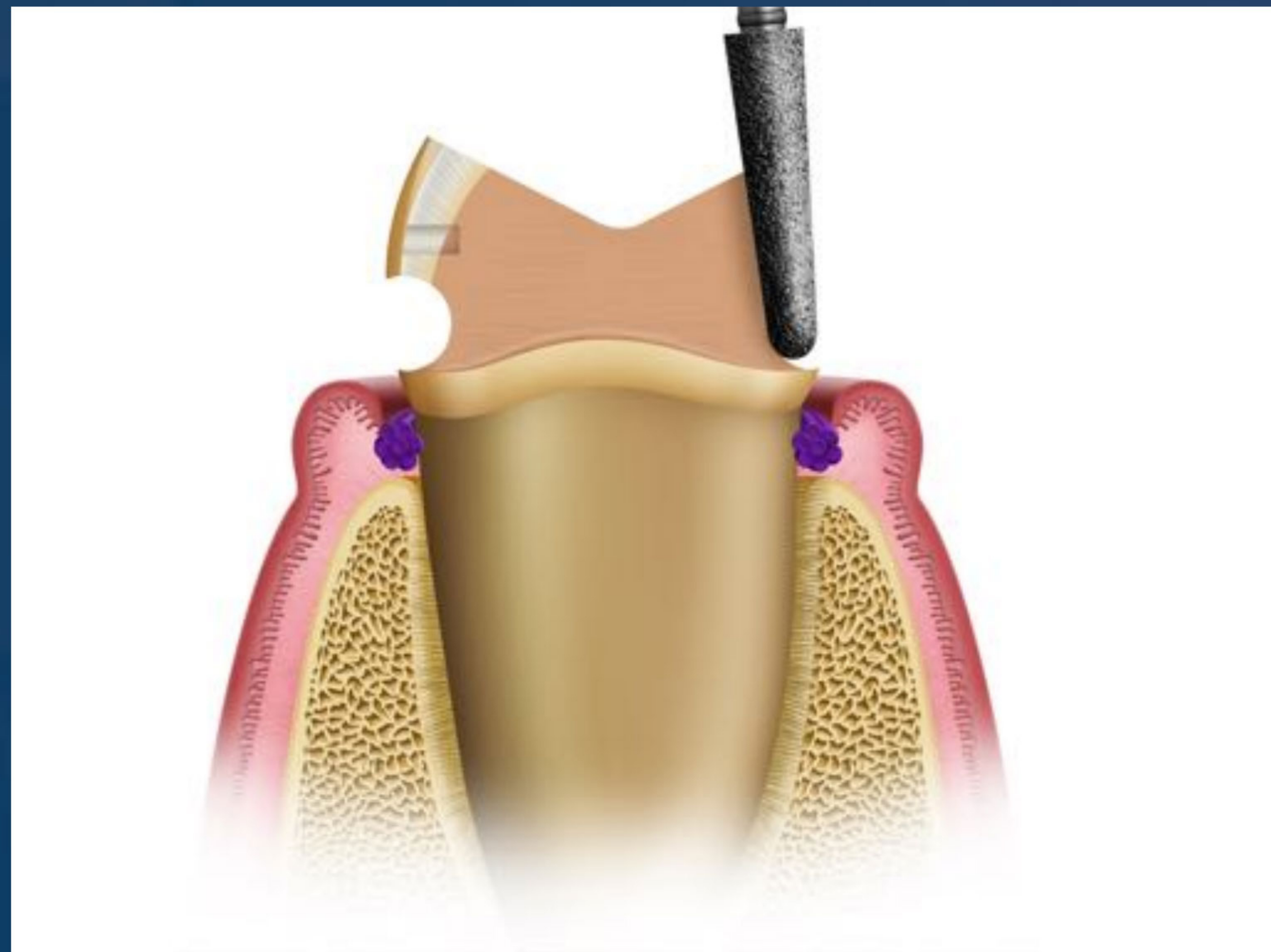


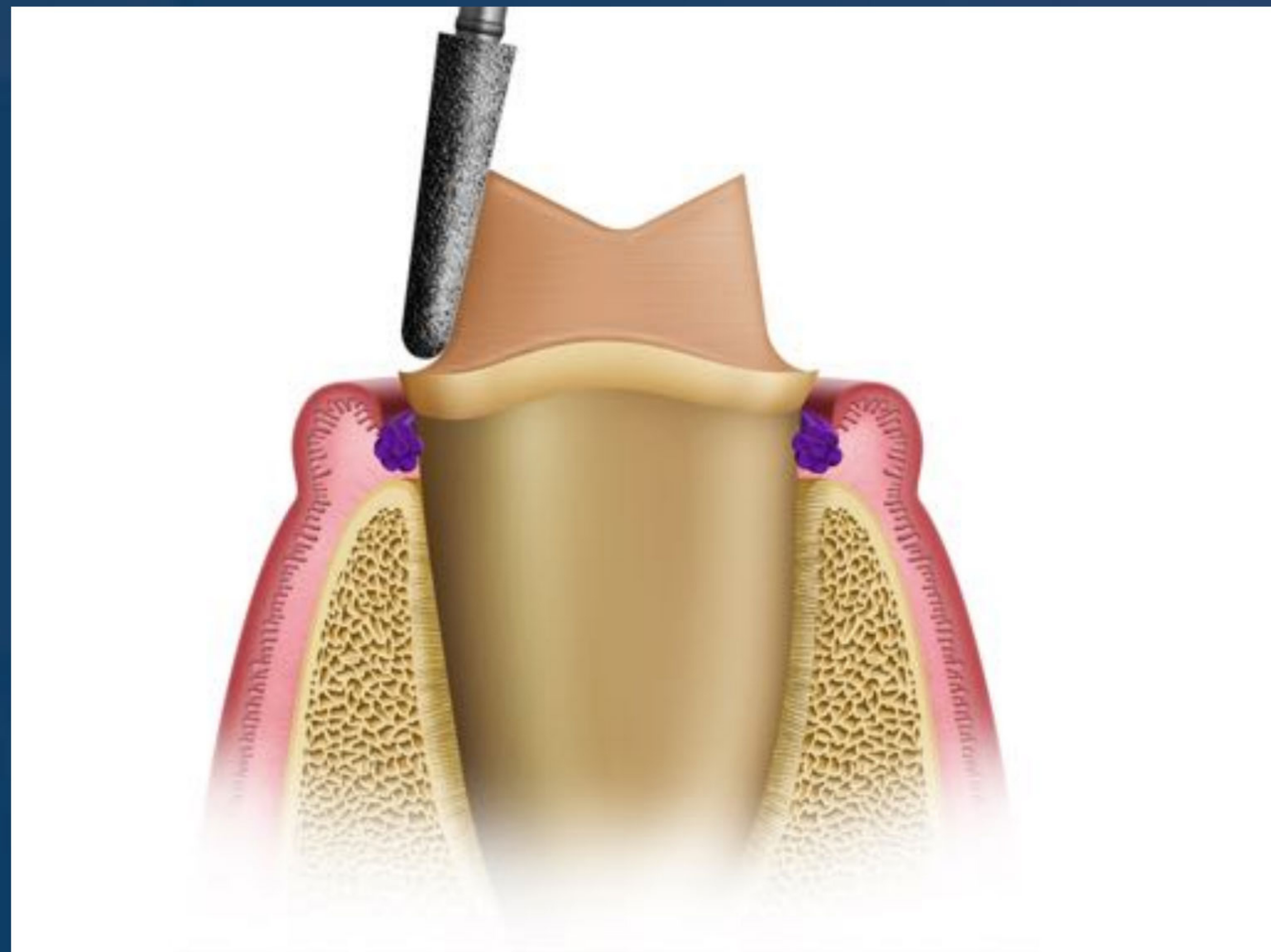


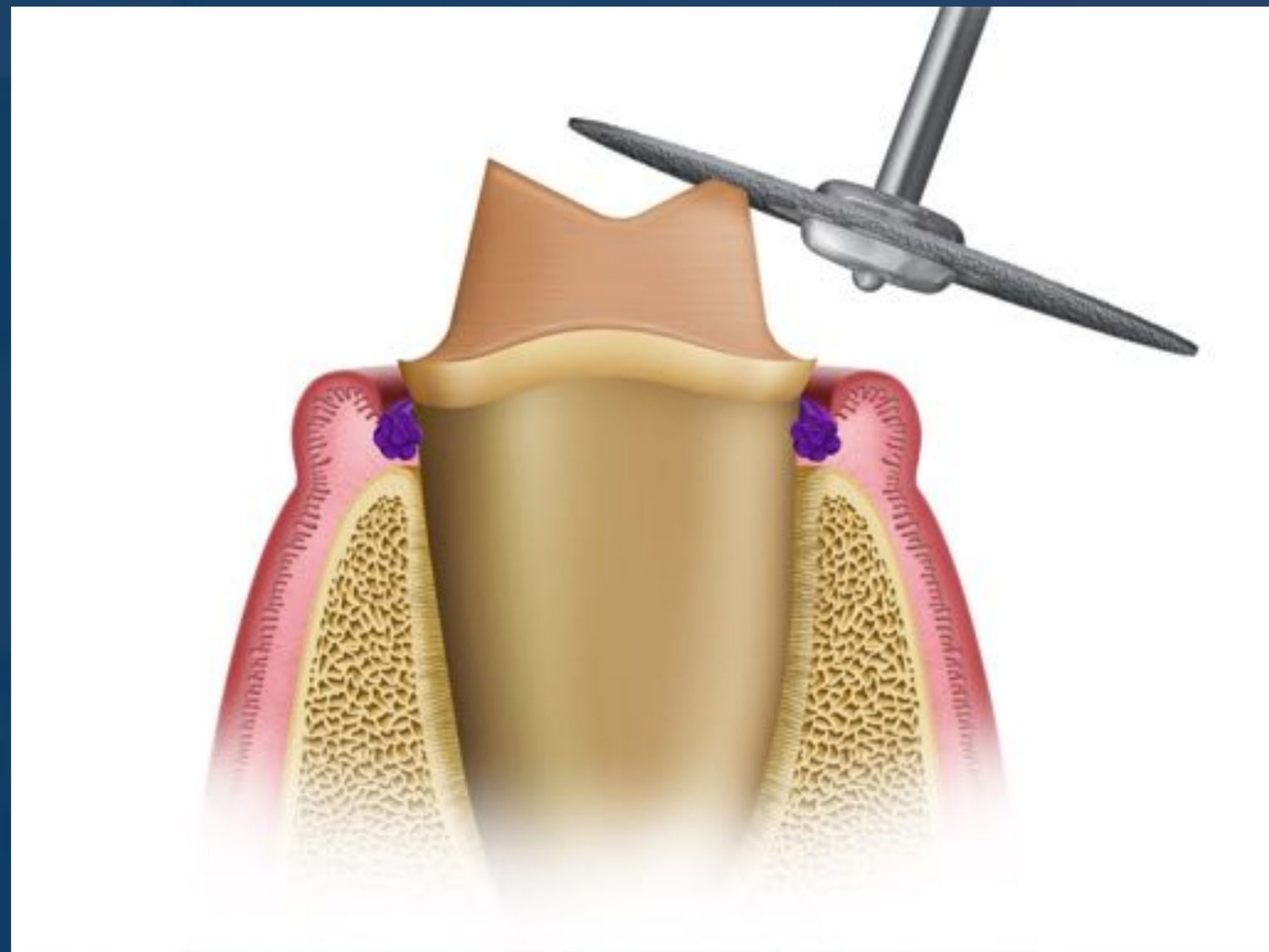

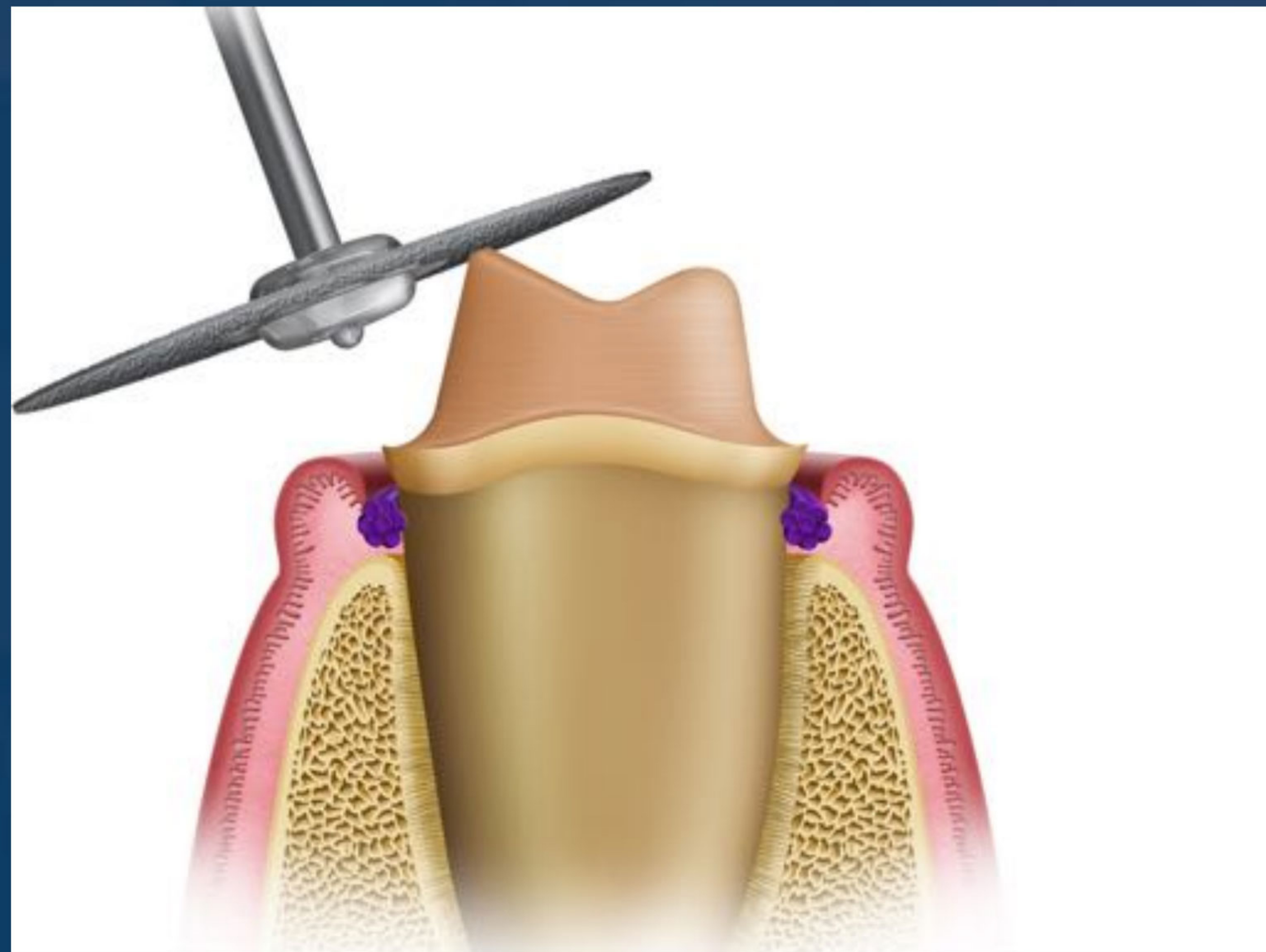



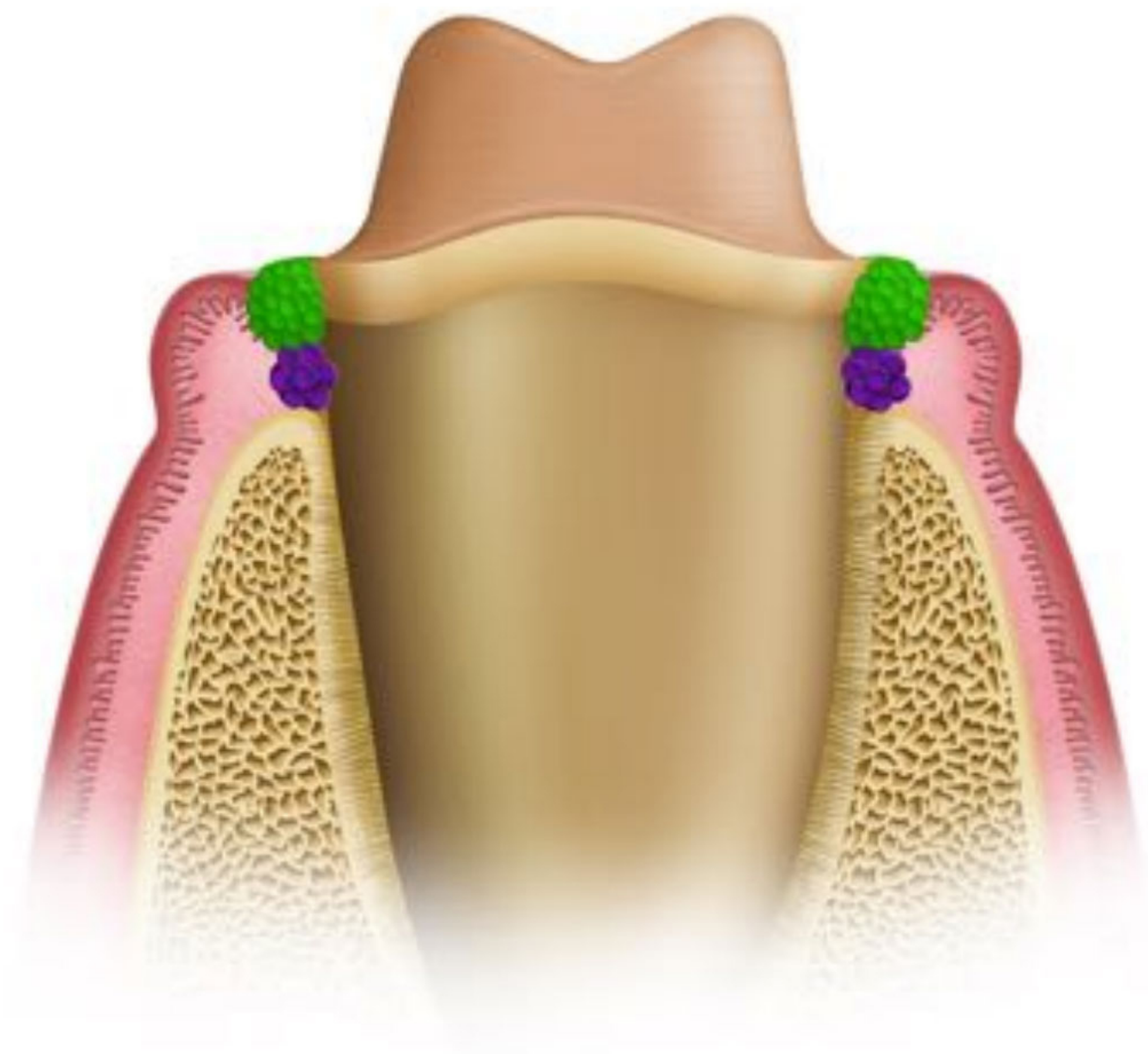

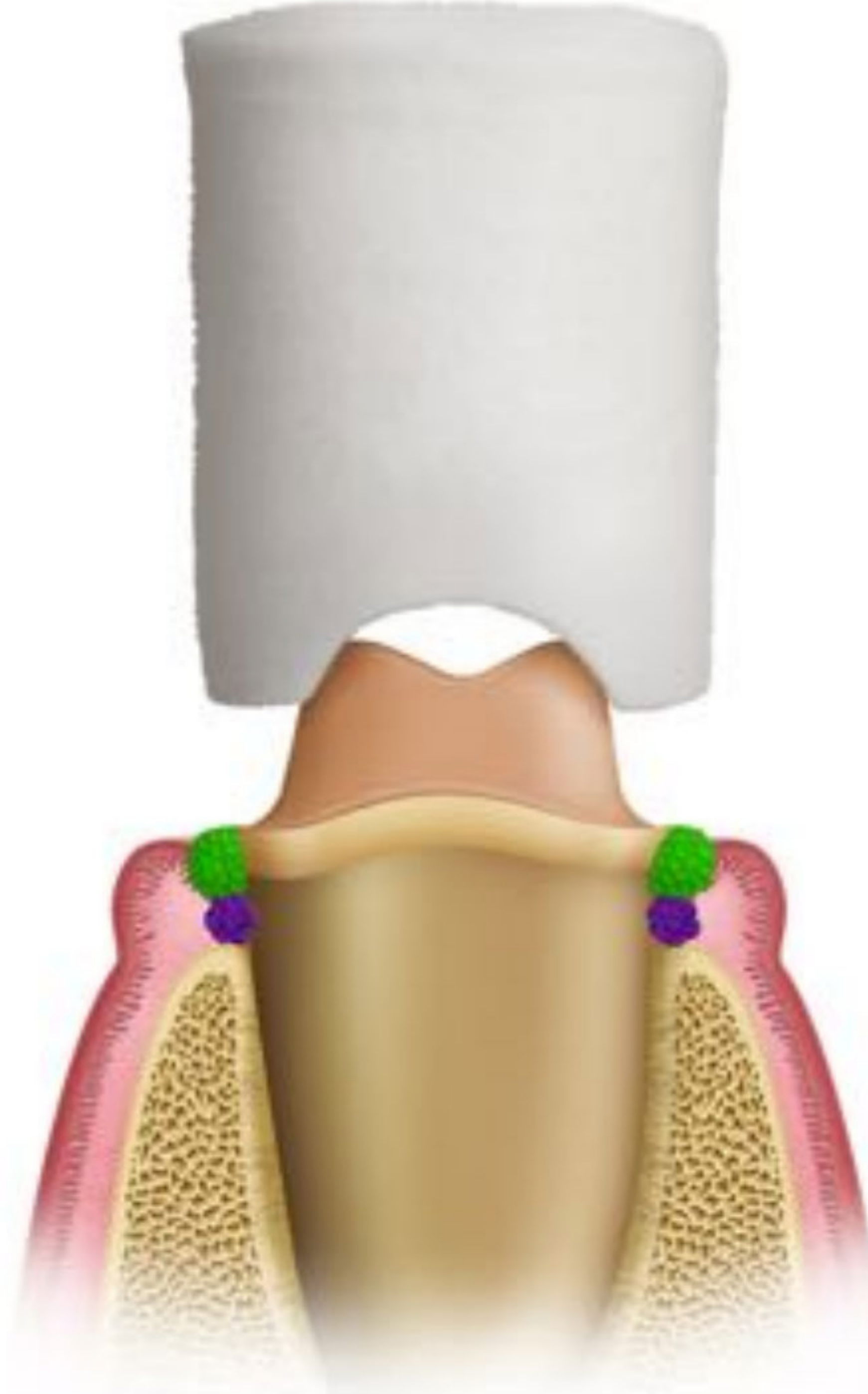



Anteriors



Step 1—Break contacts with 55 or 56 bur



Step 2—Pack plain 00 retraction cord



Step 3—Prep margin with 801 018, 016 or 014 bur



Step 4—Place incisal edge depth cuts




Step 5—Place axial reduction depth cuts



Step 6—Connect depth cuts with 856 025 SC




Step 7—Smooth margin with 856 025F

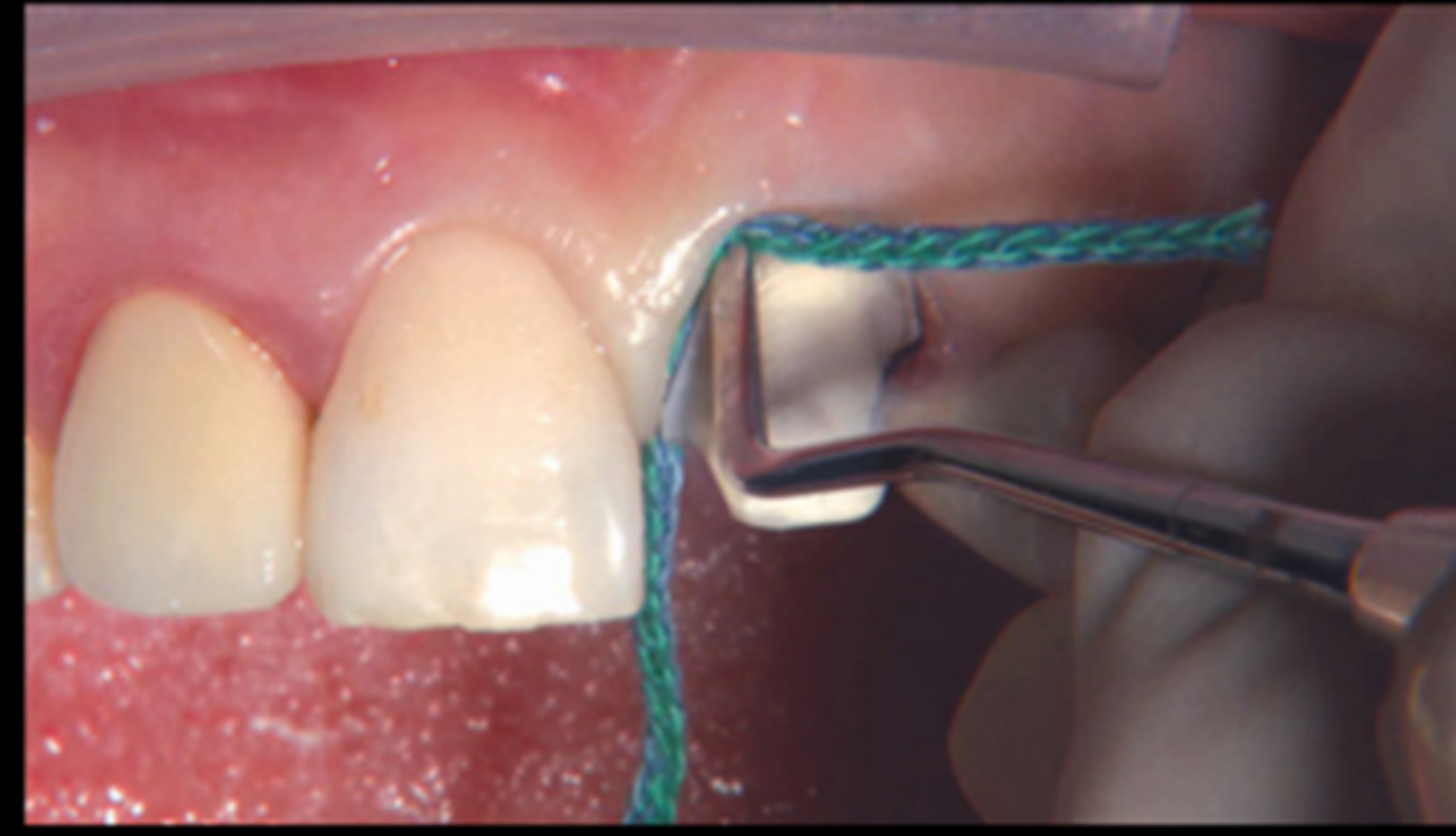

Step 8—Place top cord (size 2E or 1E)



Step 9—Place moist comprecap and bite

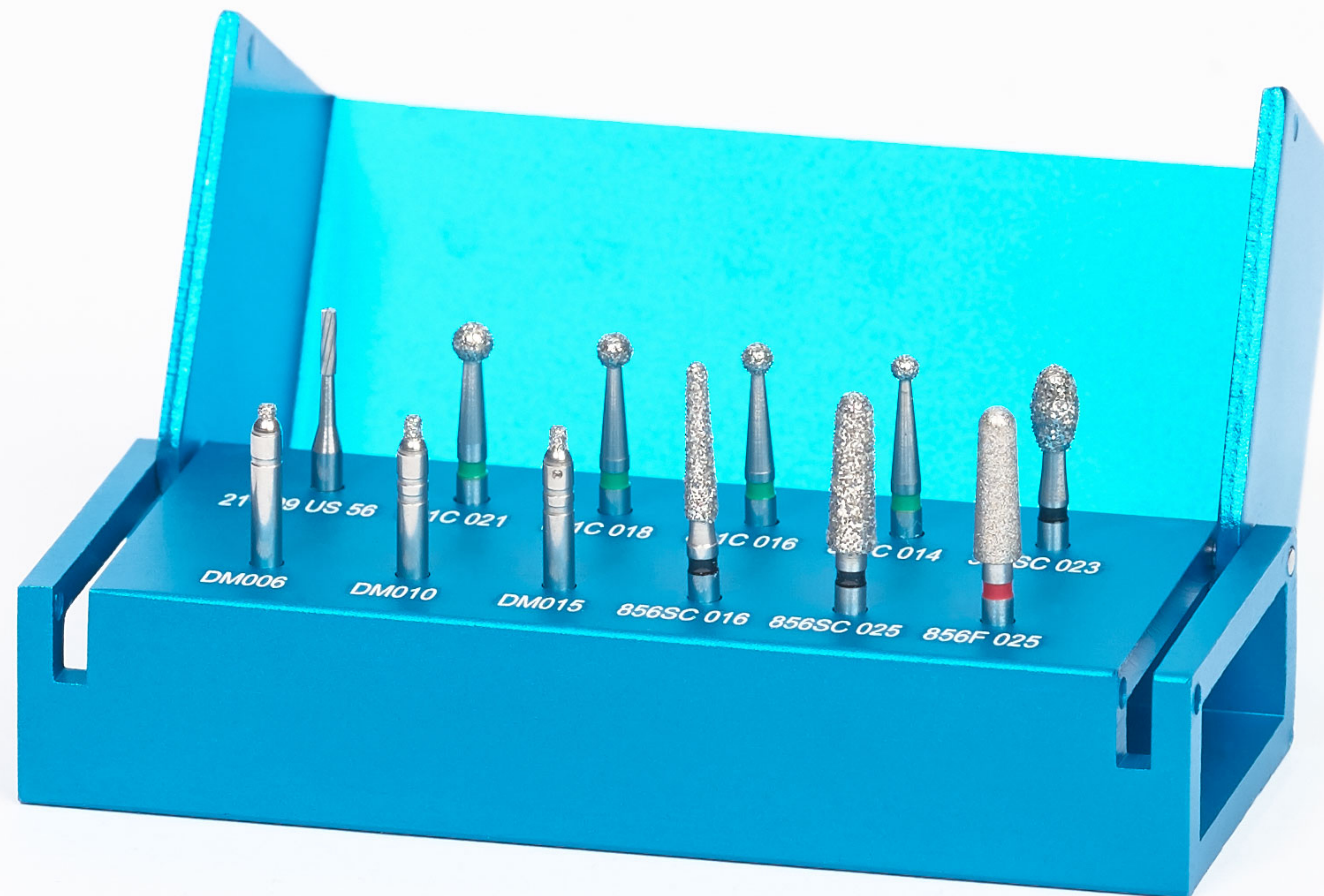

Dr. Mike's Reverse Preparation Kit

Only available through DCC

SKU: 01-001-01-DCC